




Revolutionary Design
Extraordinary Advantages



No Cannulation • 14F Delivery • Quadrafix-Seal • Repositionable

Simplified Delivery Precise Placement

The Altura™ Endograft System represents a revolutionary approach to EVAR. Its unprecedented design is comprised of a unique, repositionable bilateral D-stent system which requires no contralateral limb cannulation, delivered in a 14F ultra low-profile catheter. A broad range of patients can be treated with only six product codes, which simplifies case planning and inventory management.

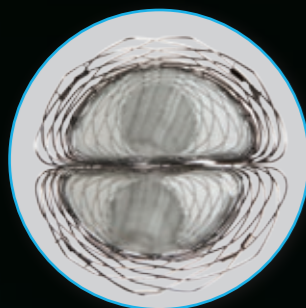


Unique Bilateral D-Stent Design

- Eliminates contralateral limb cannulation, simplifying procedures
- Braided nitinol stent provides exceptional flexibility and vessel conformity
- Ultra low-porosity woven polyester graft with endoskeleton enhances sealing properties
- Radiopaque markers for precise positioning and easy septum alignment
- Two sets of suprarenal anchors on each aortic stent prevent migration




Ultra Low-Profile 14F Delivery System

- Designed to improve access and minimise vessel trauma
- No redocking necessary for atraumatic removal
- Streamlined endograft delivery ensures controlled and accurate placement
- Built-in contrast injection ports simplify the procedure

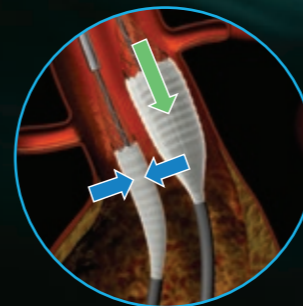


Repositionable for Accurate Placement

- Bilateral endografts may be repositioned for precise placement
- Intuitive five-step aortic endograft delivery system provides controlled deployment
- Unique one-step retrograde iliac limb deployment, for placement accuracy at the internal iliac artery

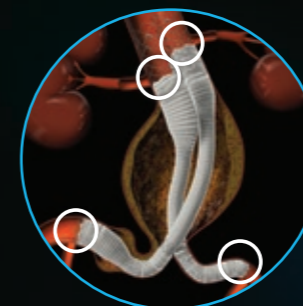


Quadrafix-Seal

- Targets each vessel and provides four points of accuracy
- Ability to independently position and offset aortic endografts maximises seal in infrarenal neck

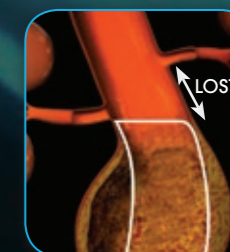

Traditional

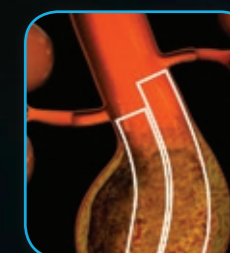

Altura

Altura's unique design eliminates gate cannulation, which can be challenging and time-consuming with a traditional stent. Coupled with the retrograde iliac deployment, these innovations make my procedures predictable, extremely precise, and easier.

– Prof. Dainis Krievins
Vascular Surgeon, Director, Institute of Science,
Stradins University Hospital, Riga, Latvia

Streamlined Case Planning Simplified Inventory Control



The Altura Endograft portfolio consists of six product codes, which allow treatment of a broad range of patients, along with simplified case planning and inventory management.

Aortic Endograft			
Product Code	Endograft Diameter (mm)	Treatment Range* (mm)	Suprarenal Stent Diameter (mm)
01-AA-24-090	24	18-22	29
01-AA-27-090	27	21-25	33
01-AA-30-090	30	24-28	36

Iliac Endograft		
Product Code	Distal Diameter (mm)	Treatment Range* (mm)
01-AI-13-065	13	8-11
01-AI-17-065	17	11-15
01-AI-21-065	21	15-18

* Treatment diameter measured from vessel lumen

Corporate Office

Lombard Medical, Inc.
6440 Oak Canyon, Suite 200
Irvine, CA 92618 | Tel: 866 790 1447
E-mail: USinfo@lombardmedical.com
LombardMedical.com

United Kingdom

Lombard Medical, Inc.
Lombard Medical House
4 Trident Park Didcot
OX11 7HJ | Tel: +44 (0) 1235 750800
E-mail: globalinfo@lombardmedical.com

Refer to product IFU for full indications, contraindications, warnings, and precautions. Altura Endograft System is CE marked and not for distribution in the USA.

© 2016 Lombard Medical Technologies, Inc. The Lombard Medical logo, the Altura Endograft System logo, and the Altura name are trademarks or registered trademarks of Lombard Medical Technologies, Inc. PN03916

